

PROJECT 11-SP-11X

N tebook No 181052


Continued From Page

TITLE: PREPARATION OF 2xMBI-11 CONTAINING A SPACER REGION

PURPOSE:

To insert a previously prepared cassette containing MBI11 (no stop codon) followed by a 8 amino acid (25 base pair) fragment, the spacer. This particular spacer is a natural sequence from *Apis mellifera* (honeybee); it occurs between identical copies of genes producing apidacins. Hopefully the spacer will have a similar effect to its natural version - to allow expression of multiple tandem copies of toxic genes. The cassette will be inserted into pET21CBDX-11 - i.e., an already-expressing vector with a copy of MBI11 that has a stop codon at the end. The cassette will be inserted such that the MBI11 with the stop codon will remain the last copy in the tandem sequence.

The final product will be as follows:

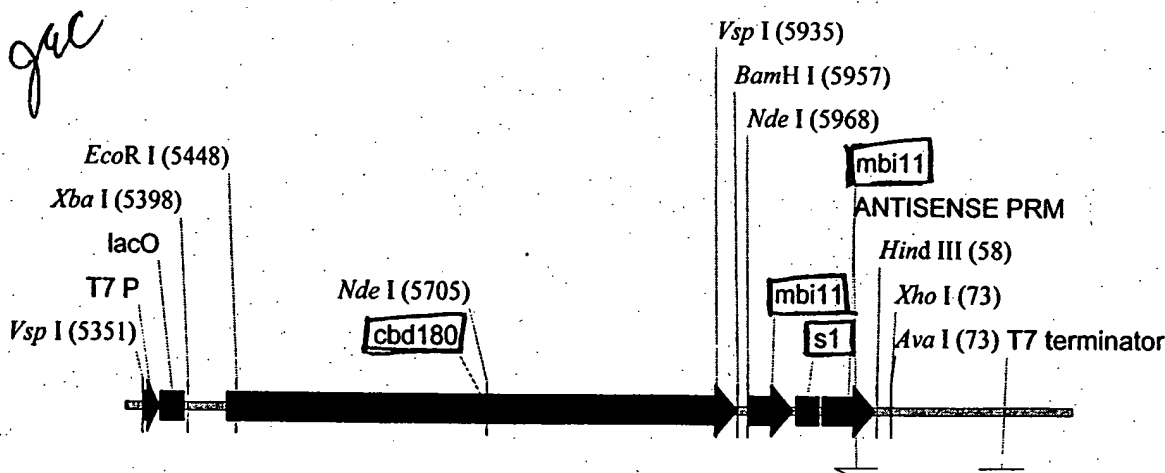

Fragment of pET21CBDX-2X11s1
965 bp (molecule 6042 bp)

Continued on Page 93

Read and Understood By

Cabralda

Signed

Date

Signed

Date

Testing expression of poly mbi-11 clones in enriched media

Strategy: Want to compare the expression of PET21CDBX-1XII w/ the other mbi-11 clones. Will see if expression decreases w/ increasing number of mbi-11 additions. It seems plausible that by enriching the media w/ lysine & tryptophan that it may be possible to increase expression, especially, the mbi-11 clones. (Since the amino acid sequence of mbi-11 is rich in lysine & tryptophan). Will ^{compare} expression by optimizing the amount that is loaded on the SDS-PAGE gels.

Methods:

25mg/ml L-lysine stock solution: 250mg L-lysine (lot #35H12175) dissolved in 10ml sterile distilled water.
- filter sterilized thru 0.2um filter

12.5mg/ml L-tryptophan stock solution: 250mg L-tryptophan (lot #73H035225) (con 0.1M HCl) in 19.8ml distilled, 172uM conc. HCl.
- filter sterilized thru 0.2um filter

Enriched Tergitol broth - add 18uM of stock lysine soln & 36uM of stock tryptophan soln to 18uM TB.

Inoculate 3ml TB (w/ 100ug/ml amp + 30ug/ml Kan) w/ PET21CDBX-1XII, or PET21CDBX-2XII, or PET21CDBX-3XII, or PET21CDBX-4XII. (w/ in XL-1 + pGP1-2).
Incubate @ 30°C, w/ shaking.

Also inoculate 2ml TB w/ 300ug/ml Kanamycin w/ XL-1 + pGP1-2.

Prepare 4x 3ml enriched TB w/ 100ug/ml amp + 30ug/ml Kan.
Inoculate w/ PET21CDBX-1XII/2XII/3XII/4XII in XL-1 + pGP1-2.

Also inoculate 2ml enriched TB w/ 300ug/ml Kan w/ XL-1 + pGP1-2.
Incubate @ 30°C, w/ shaking.

Continued on Page 83

Read and Understood By

Spulman
Signed

Date

Di-
Signed

Date

Induce expression of λ phage overnight culture grown @ 30°C by adding 1ml fresh TB (w/ and w/out lysine + tryptophan); shaking @ 42°C, 31 hours.

Test Abcd of overnight cultures: (16x dilution)

Enriched media:

XL-1 + pGP1-2 : 0.603

PETZICBDX-1X11 : 0.445

" - 2X11 : 0.091

" - 3X11 : 0.339

" - 4X11 : 0.081

Normal media:

Enriched media: XL-1 + pGP1-2 : 0.541

PETZICBDX-1X11 : 0.476

" - 2X11 : 0.111

" - 3X11 : 0.345

" - 4X11 : 0.093

Test Abcd of induced cultures: (16x dilution)

Normal media:

~~XL-1 + pGP1-2~~ PETZICBDX-1X11 : 0.430

" - 2X11 : 0.091

" - 3X11 : 0.231

" - 4X11 : 0.138

Enriched media: PETZICBDX-1X11 : 0.425

" - 2X11 : 0.139

" - 3X11 : 0.220

" - 4X11 : 0.103

* Spin down cells equivalent to OD = 1.5. Resuspend in 100ul SDS-Reducing buffer. Boil samples, 10'. Load 10ul = O.D. = 0.15 of each sample onto 15% SDS-PAGE.

Note: PETZICBDX-4X11 - 42°C sample, must have accumulated under, during boiling step. Loaded 20ul (eyeball estimate). - prepared another sample for use next time.

* Abcd, enriched to equal OD = 1.5 (then resuspended in 100ul reducing buffer)

Overnight cultures:

Enriched media:

Induced cultures:

Enriched media:

XL-1 + pGP1-2 - 173ul

XL-1 + pGP1-2 - 155ul

PETZICBDX-1X11 - 218ul

1X11 - 221ul

PETZICBDX-1X11 - 197ul

1X11 - 209ul

2X11 - 966ul

2X11 - 674ul

" 2X11 - 845ul

2X11 - 1.030ml

3X11 - 466ul

3X11 - 426ul

" 3X11 - 272ul

3X11 - 271ul

4X11 - 679ul

4X11 - 910ul

" 4X11 - 1.008ml

4X11 - 1.157ml

Continued on Page 84

Read and Understood By



Haleman
Signed

Date

[Signature]
Signed

Date


Poly mbi-

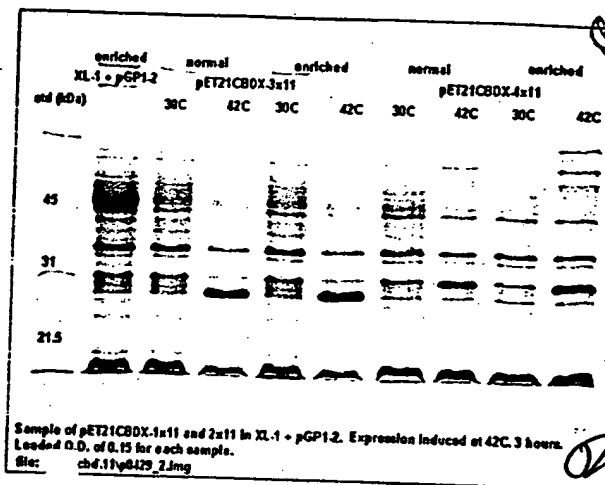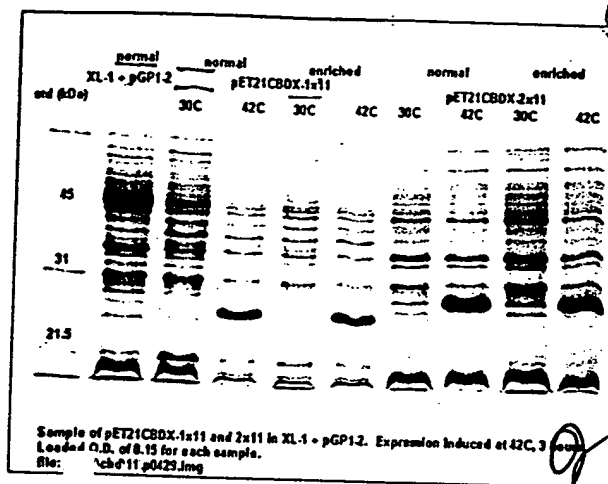

Good expression.

No protein analysis on protein gels.

pET21CBDX-1x11:

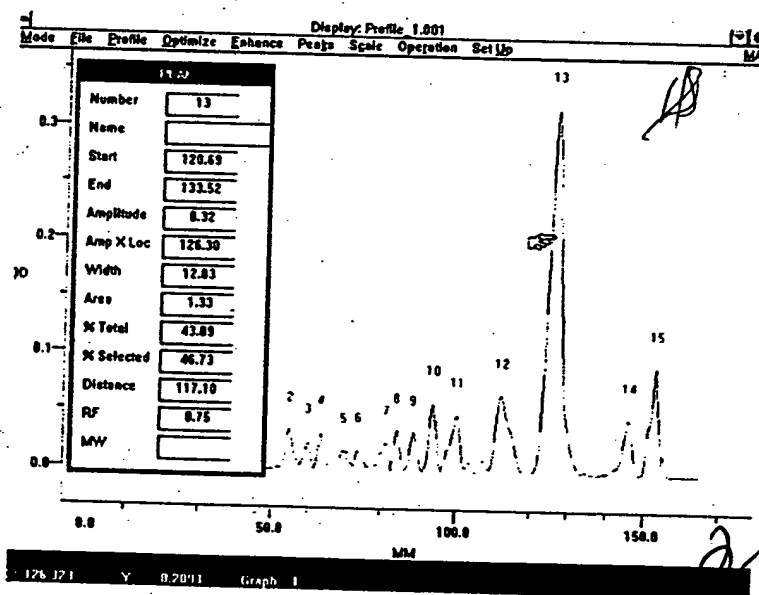

CBDX = 180 amino acids
1x11 = 13 aa.

*fusion protein = 43%
of total protein*

*1x11 = 6.7% peptide in
fusion protein*

Continued on Page

85

Read and Understood By


[Signature]
Signed

Date

[Signature]
Signed

Date


PETRICBDX-ZXII

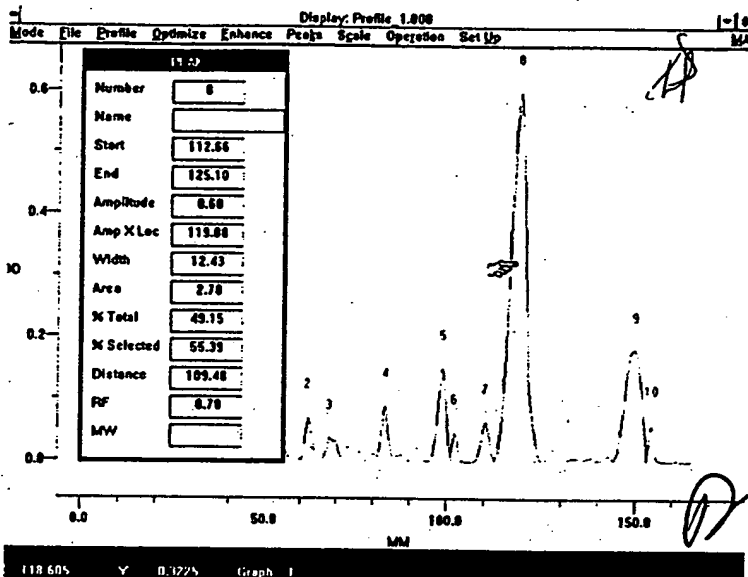

fusion protein = 49% of total protein

ZXII = 13.6% peptide in fusion protein

PETRICBDX-ZXII

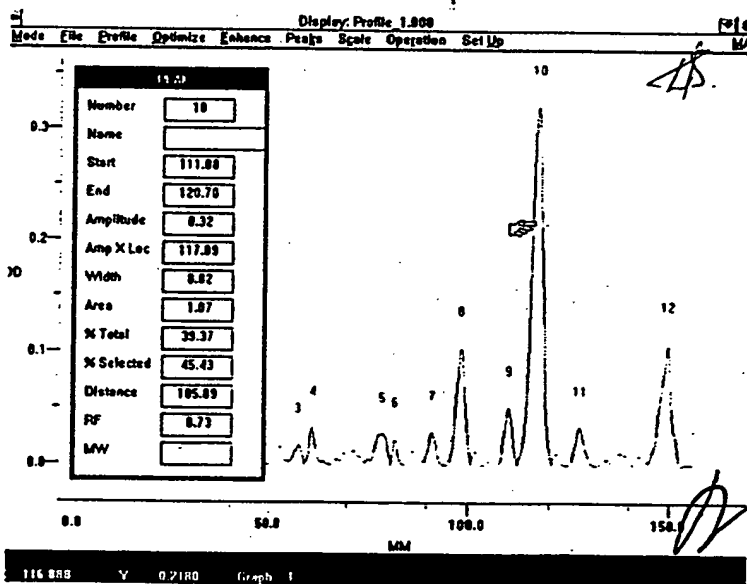

fusion protein = 39% of total protein

ZXII = 13.8% of fusion protein = peptide

Continued on Page 86

Read and Understood By


Signed

Date

Signed

Date

PET21CDBX-4x11

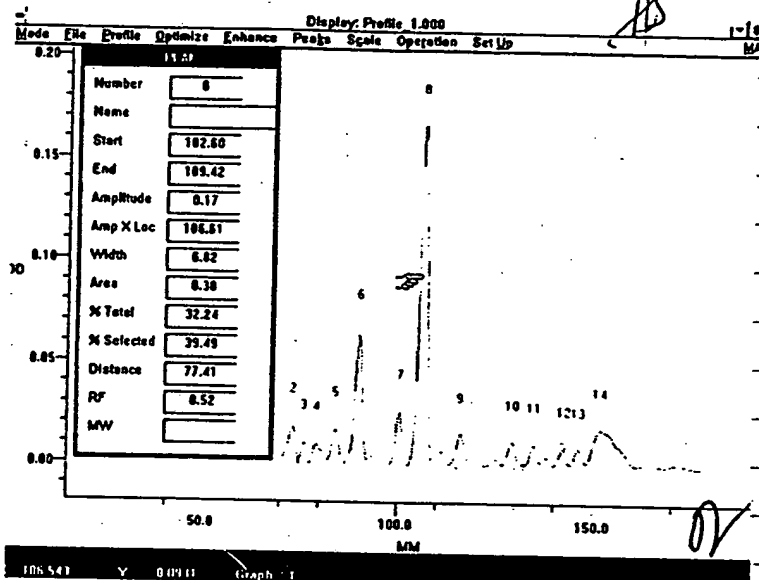

fusion protein = 32% of
total protein

peptide = 22.4% of fusion
protein

But same band was found
@ 30°C (uninduced sample)

→ will redo expression of
PET21CDBX-4x11.

Summary:

| mbi-11 clone | % expression of fusion protein | % peptide | % peptide expression |
|--------------|-----------------------------------|-------------|-------------------------|
| 1x11 | 43.09 | 6.735751295 | 2.902435233 |
| 2x11 | 49.15 | 12.62135922 | 6.203398058 |
| 3x11 | 39.37 | 17.80821918 | 7.01109589 |
| 4x11 | 32.24* | 22.4137931 | 7.226206897* |

* percentages are lower when the 30°C (uninduced) samples
are taken into consideration

Continued on Page

Read and Understood By

Signed

Date

Signed

Date

TITLE: EXPRESSION OF CBDX-2X11S1 FUSION PROTEIN IN E. COLI
MC4100(F) AND BL21(DE3)

PURPOSE:

Two strains were previously prepared:

- 1.) MC4100(F) + pGP1-2 + pET21CBDX-2X11S1
 - 2.) BL21(DE3) + pET21CBDX-2X11S1
- (see pp. 92-100, MB105A, and pp. 1-22, MB1082).

Both of these strains should be capable of expressing the CBDX-2X11S1 fusion protein via the T7 promoter system (in MC4100(F), after induction at 42°C; in BL21, constitutively with overnight growth).

In this experiment I will test fusion protein expression from each of these strains using the appropriate induction methods.

EXPERIMENTAL DETAILS:

I GROWTH OF STRAINS

Overnight Growth

- four culture tubes were prepared: 2 containing 2 mL TB + 100 µg/mL Amp, and 2 containing 2 mL TB + 100 µg/mL Amp + 30 µg/mL Kan.
- the Amp tubes were both inoculated with the plasmid prep O/N culture (see p. 22) of sample #2 (BL21)
- the Amp + Kan tubes were both inoculated with the plasmid prep O/N culture (see p. 22) of sample #11 (MC4100(F))
- tubes were placed @ 37°C / 30°C O/N with shaking


Continued on Page 28

Read and Understood By

J. Cabralda

Signed

Date

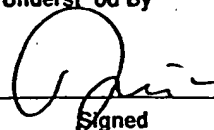

Signed

Date

- note: these O/N cultures from plasmid preps are quite old (several days), so I'm preparing 2 cultures of each and will test for presence of plasmids following O/N growth.

HARVESTING OF BL21 SAMPLES

Harvesting

- removed 100 μ l from BL21 O/N's \rightarrow & 13,200 rpm 1 min
- \rightarrow removed S/N
- resuspended pellet in 100 μ l SDS buffer
- froze @ -20°C .

INDUCTION OF MC4100 (F) SAMPLES

Removal of Preinduction Samples

- removed 100 μ l from MC4100 (F) O/N's and harvested as described above \rightarrow these are the preinduction samples
- resuspended in 100 μ l + froze @ -20°C .

Induction

- removed 0.9 ml of each MC4100 (F) O/N into a fresh sterile culture tube
- added 0.9 ml TB to each
- placed @ 42°C for 3h.
- following induction, removed 200 μ l \rightarrow harvested as described above
- resuspended in 100 μ l \rightarrow froze @ -20°C .

Continued on Page 29

Read and Understood By

J. Cabralda

Signed

Date

Jm

Signed


Date

IV PLASMID PREPSPlasmid Preps

- during induction time, removed remainder of all 4 O/N cultures to microcentrifuge tubes (~1 mL) and performed plasmid preps as per SOP.
- stopped after NH_4OAc precipitation
- resuspended in 50 μL TE + RNase.

Gel

- ran 1 μL of each prep with 1 μL pGPI-2 + pET21CDBX-2X11S1
- ran on 1% agarose for 45' @ 80V.

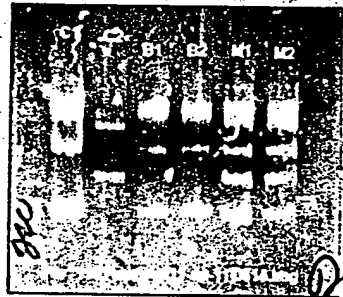

C1 = pET21CDBX-2X11S1
C2 = pGPI-2.

- appropriate plasmids are all there!
- need only test expression of one of each pair

V TESTING OF EXPRESSION LEVELPreparation of Samples

- boiled samples for 10 min.
- spun down

Gel

- prepared ⁴⁰⁰ 15% / 4% stacking gel
- loaded 10 μL of BL21 (1) and Mc4100 (K) (1), along with an expressing 2x sample from Jan

Continued on Page 30

Read and Understood By


J. Cabralda
Signed

Dat

J. Jan
Signed

Dat



run @ 100V for ~ 100 min

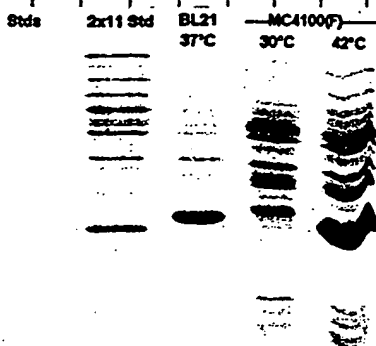

gac

File: lapacer11-ap-110616.img

Profile Analysis:

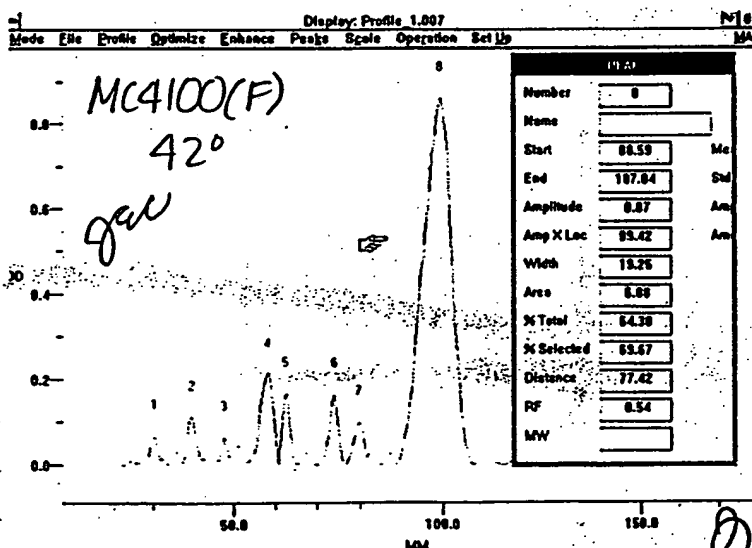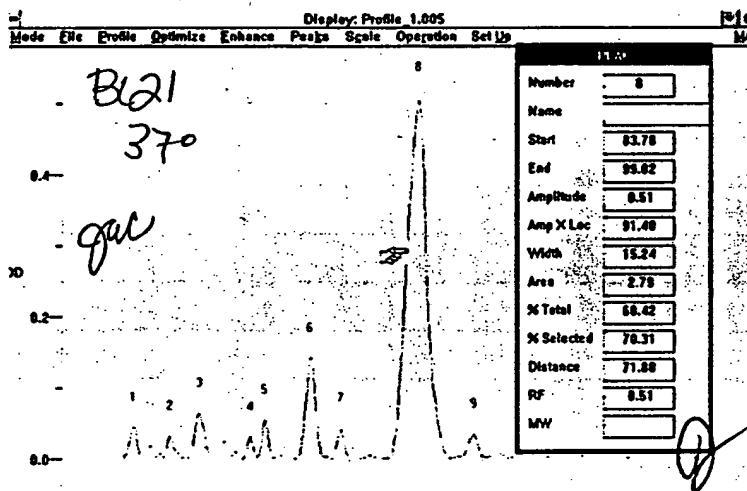

Continued on Page 31

Read and Understood by

J. Cabralda
Signed

Date

[Signature]
Signed

Date

PROJECT 2X11S1 EXPRESSION

Continued From Page

30

EXHIBIT page 11 of 23

- apparently very high expression!
- however, band appears too big - maybe a frame-shift mutation somehow?
- will do sequencing to check.
- NOTE: insert also looked too big - see p.7, fig. 9.

DISCUSSION:


Expression of 2X11S1 fusion protein, at apparently high levels, was achieved in both BL21(DE3) and MCH100(F). However, the fusion protein band appears perhaps slightly larger than expected. The clone will be sequenced.

Continued on Page

Read and Understood By

J. Cabialda
Signed

Date

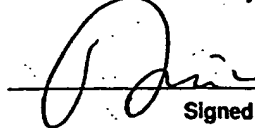
Signed

Date

TITLE: TESTING THE EXPRESSION OF CBD96-4x11 AND CBD96-2x11B7
FUSION PROTEIN IN E. COLI MC4100(F)

PURPOSE:

In a previous experiment (see pages 35-44, this book) the expression of CBD96-6x11 and -10x11 fusion protein was tested and determined to be nil. In this experiment, CBD96-4x11 will be tested for expression to see whether it expresses at a greater level than CBDX-4x11.

In addition, Jan has prepared a CBD96-2x11B7 clone. This clone was prepared from a 2x11B7 fragment that was prepared in one piece - i.e., it wasn't constructed from two 1x clones. MB11B7 is, as of now, our "official lead compound" so expression at high levels is essential. The expression of this fusion protein will also be tested.

EXPERIMENTAL DETAILS:

I CREATION OF TEST STRAINS

Electroporation

- used MC4100 (F) as the host strain
- chose clones #8 + 12 of 2x11B7, and #4-4 and 4-6 of 4x11
- electroporated plasmid DNA, from Jan Burman 1 μ l 2x11B7 clones or 2 μ l 4x11 clones (low [DNA]) with 1 μ l p6PI-2.
- plated 10 μ l and 50 μ l on Mac + Amp + Kan
- grew O/N at 30°C.
- few colonies - especially for 4x11 clones
- but at least one colony on each 50 μ l plate

Continued on Page 45

Read and Understood By

J. Cabralda

Signed

Date

J. Burman

Signed

Date

II OVERNIGHT GROWTH OF STRAINSGrowth

- prepared four culture tubes containing 3 mL TB + 100 µg/mL ampicillin + 30 µg/mL kanamycin
- inoculated each tube from with a single colony from one of the 50 µL electroporation plates
- grew O/N at 30°C with shaking (300 rpm).

III PREINDUCTION HARVESTING & INDUCTION OF EXPRESSIONPreinduction Harvesting

- forgot to remove preinduction samples!
- so added 100 µL fresh TB to tubes after they were empty and so grew at 30°C for ~3h.
- then harvested (100 µL) → & → removed S/N
- resuspended in 100 µL SDS reducing buffer.

Induction of Expression

- removed 1 mL of each sample and put in a tube containing fresh TB
- incubated @ 42°C with shaking (300 rpm) for 3h.

IV CHECKING FOR PRESENCE OF PLASMIDPlasmid Preps

- meanwhile, while samples are expressing, use rest of O/N cultures (2 mL) to do plasmid preps
- did as per SOP to NH_4OAc precipitation
- resuspended pellets in 50 µL TE + RNase

Continued on Page

43

Read and Understood By

J. Cabralda

Signed

Date

[Signature]


Signed

Date

Gel

- prepared a 1% agarose gel
- loaded 3 μ l pET21CB096 (from GPN) and 3 μ l pGP1-2 as controls
- loaded 1 μ l each plasmid prep and 5 μ l 1 kb ladder.
- ran 45' @ 80V.

Figure 1 Agarose gel of putative transformants. All clones are positive for the appropriate plasmids.

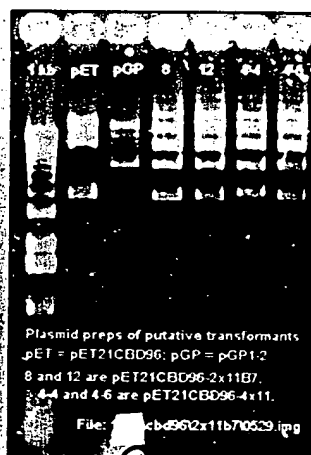

V CHECKING FUSION PROTEIN EXPRESSION

Harvesting of Expressed Plasm Samples

- removed samples from 42°C incubation after 3h.
- extracted 200 μ l per tube \rightarrow ϕ \rightarrow removed GPN
- resuspended each pellet in 100 μ l SDS buffer.

Preparation of Samples

- removed all samples (30°C and 42°C) from storage at -20°C.

Continued on Page

49

Read and Understood By

J. Cabralda

Signed

Dat

J. Cabralda

Signed


Dat

- boiled samples for 10 min, then spun down.

Gel

- prepared a 15%/4% stacking PA gel
- loaded 10 μ l of each sample with 5 μ l low molecular weight protein standards.
- ran 45' @ 200V.

Figure 2 Expression of 2x11B7 and 4x11 from CBD96.

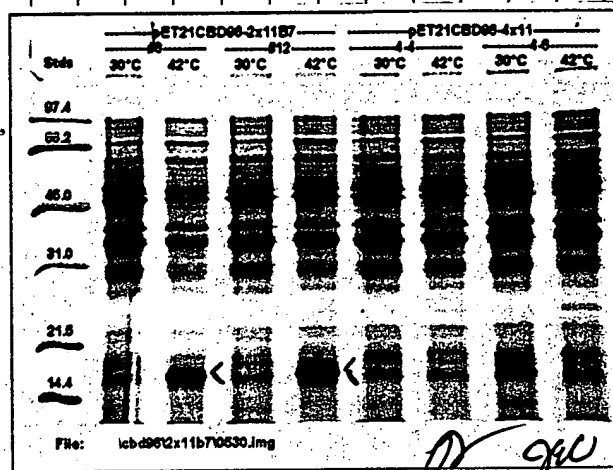

- since 2x11B7 is being expressed at fairly low levels and 4x11 isn't expressed at all, I'll check to see whether any/all of the fusion protein is soluble in the cell as opposed to sequestered in inclusion bodies.

VI OVERNIGHT GROWTH

Growth

- set up four new culture tubes containing 25mL TB + 100 μ g/mL ampicillin + 30 μ g/mL kanamycin
- inoculated using loopfuls of last night's overnights (30°C).

Continued on Page 51

Read and Understood By

G. Cabralda

Signed

Date

[Signature]

Signed

Date

- grew O/N at 30°C, with shaking (300 rpm)
- NOTE: in morning, incubator was at $\geq 32^{\circ}\text{C}$ and "Hi-Temp" alarm was flashing.

VI PREINDUCTION HARVESTING & INDUCTION OF EXPRESSION

Preinduction Harvesting

- in morning, removed ~~20~~ 100 μL of each culture to an Eppendorf tube \rightarrow ϕ \rightarrow removed S/N.
- resuspended pellet in 100 μL SDS reducing buffer
- freeze @ -20°C until needed.

Induction of Expression

- prepared 4 fresh culture tubes, each containing 0.9 mL TB
- added 0.9 mL culture to tubes
- incubated @ 42°C for 3h with shaking (300 rpm)

VII CHECKING FUSION PROTEIN EXPRESSION - SOLUBLE/INSOLUBLE

Harvesting of Postinduction Samples

- removed 200 μL of each sample after incubation
- ϕ \rightarrow removed S/N \rightarrow resuspended in 100 μL SDS reducing buffer

Sonication to Release Soluble Cell Components

- transferred 1 mL induction culture to an Eppendorf \rightarrow ϕ and removed S/N
- resuspended in 980 μL H₂O + 20 μL 0.5 M EDTA (final [EDTA] = 0.010 M) \rightarrow vortexed to mix
- Sonicated @ 10 (amplitude) using microtip: 1 minute \times 3 times (kept on ice between)


Continued on Page 52

Read and Understood By

J. Cabralda

Signed

Date

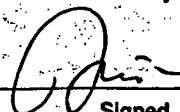

Signed

Date


Harvesting of Sonicated Samples

- removed 100 μ l each sample \rightarrow diluted w/ 100 μ l SDS reducing buffer
 \rightarrow this is the SONICATED sample
- removed & remaining sample
- resuspended pellet in 400 μ l buffer
 \rightarrow this is the INSOLUBLE sample
- transferred 200 μ l SIN to separate tube and added 200 μ l buffer
 \rightarrow this is the SOLUBLE sample

Gels

- prepared two 15%/4% stacking PA gels

Figure 3 PA gels of expression samples. 10 μ l per well.

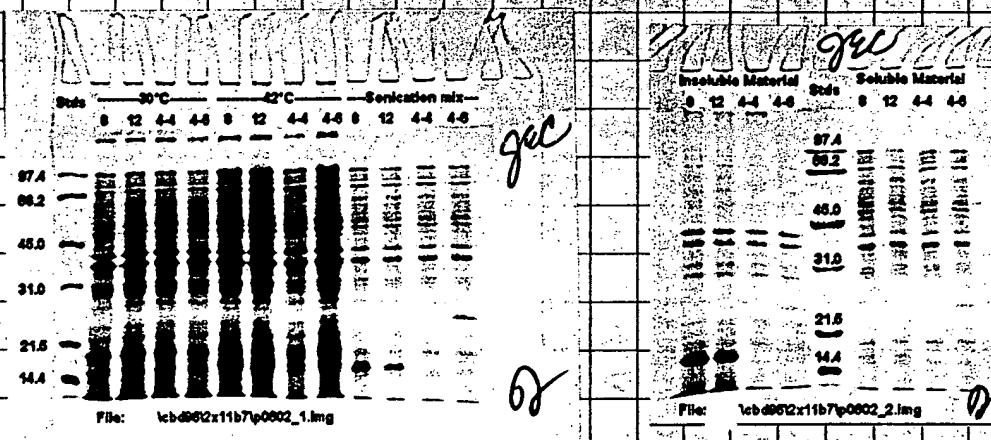

- prepared an additional 15%/4% gel and ran the 4-6 samples only, trying to "even out" loading to get ~equal amount of material per lane

Continued on Page 57

Read and Understood By

J. Cabralda
Signed

Dat


Signed

Dat

- i.e., loaded 7.5 μ l 30°C
 5 μ l 42°C
 12.5 μ l sonic.
 12.5 μ l insol.
 12.5 μ l sol.

- ran @ 200V, 45'

Figure 4 Expression samples of pET21CBD96-4x11 clone #4-b.

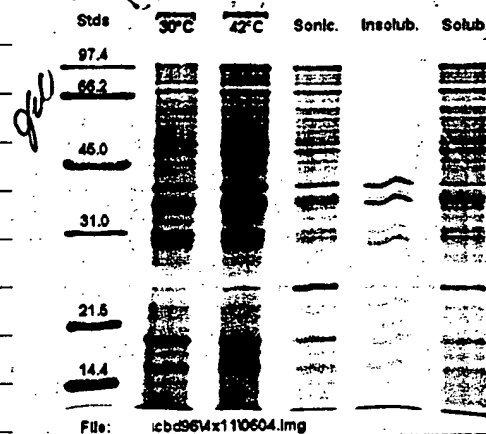

DISCUSSION:

Good expression of pET21CBD96-2x11B7 was obtained. It appears from figure 3 that CBD96-2x11B7 fusion protein is insoluble.

Interestingly, a low level of CBD96-4x11 fusion protein was expressed, but it appears to be soluble, with only a very light band appearing in the "insoluble" lane. For unknown reasons, the fusion protein band is not visible in the 42°C lane. Perhaps this sample was not completely lysed prior to loading.

If CBD96 fusion proteins are soluble, this could cause many problems in later purification steps.

C ntinu d n Pag

R ad and Understo d By

J. Cabralda

Signed

Dat

[Signature]

Signed


Dat

TITLE: TESTING EXPRESSION OF CBD96-11B7-sp-2x11B7 AND
CBD96-2x11B7-sp-2x11B7 FUSION PROTEINS

PURPOSE / BACKGROUND INFO:

I have obtained pET21CBD96-11B7-sp-2x11B7 and
pET21CBD96-2x11B7-sp-2x11B7 clones from Jan. These
are as shown in figure 1.

Figure 1 3x11B7sp and 4x11B7sp clones used in this experiment.

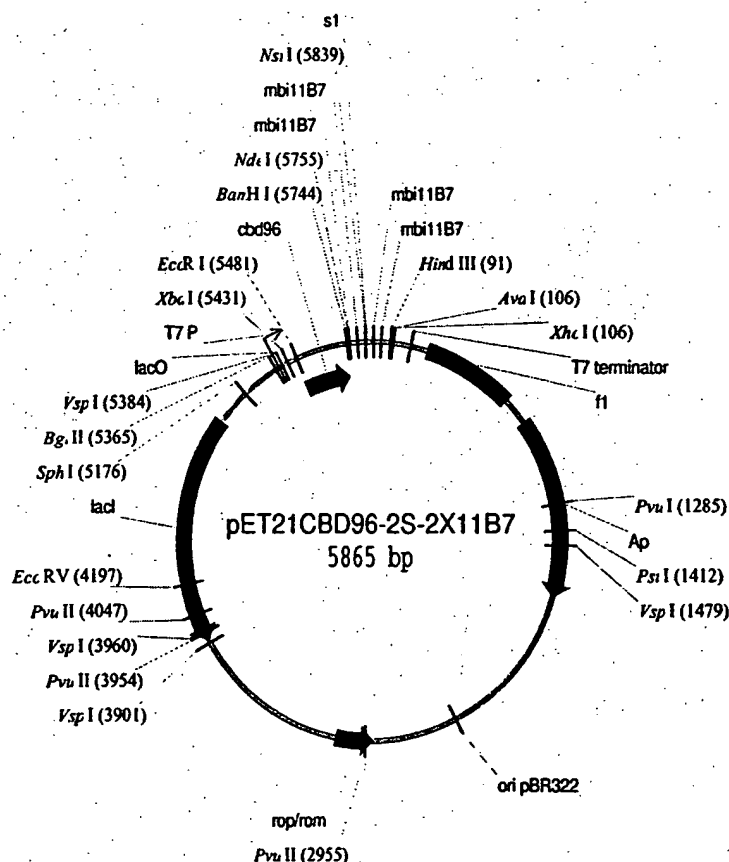

From pET21CBD96-2x11B7(LAST)
+ pBCKS-V-2x11B7S



Fragment of pET21CBD96-2S-2X11B7
474 bp (molecule 5865 bp)

Jac

OR

Continued on Page 69

Read and Understood By

J. Cabralda
Signed



Date

[Signature]
Signed

Date

From pET21CBD96-2x11B7 (AST)
+ pBCKS-V-11B7S.

gac

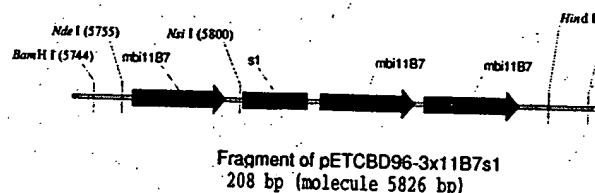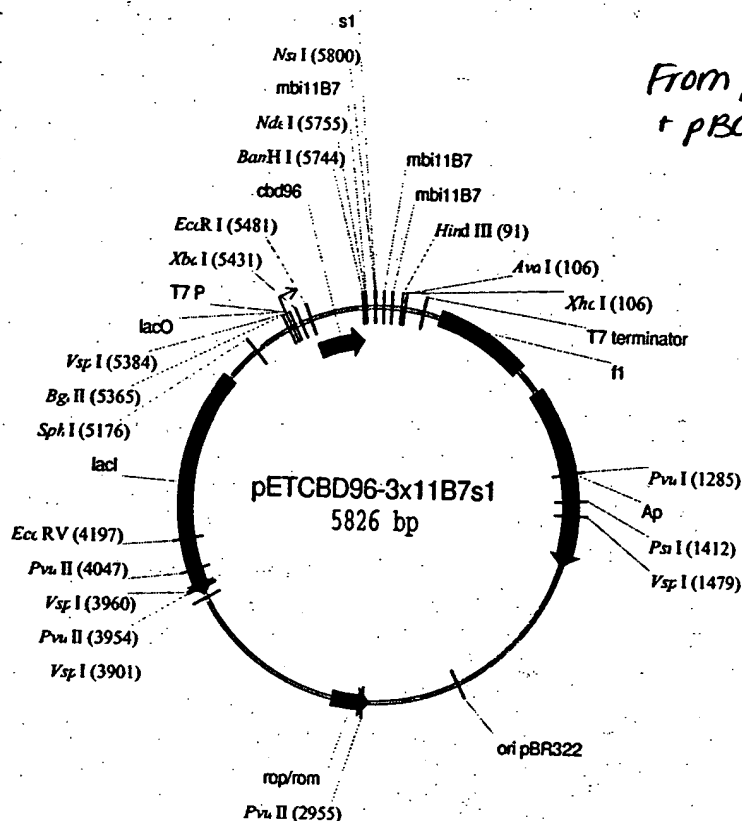

I will now create expression strains containing these clones and test for ability to express fusion protein.


EXPERIMENTAL DETAILS:

I CREATION OF EXPRESSION STRAINS

Electroporation

- obtained DNA from Jan

Figure 2 B/H digests of Jan's clones which showed inserts of correct size (i.e., all of these clones are positive).

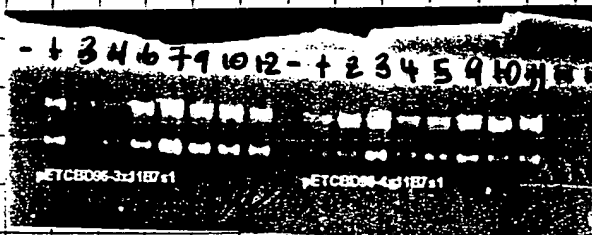

Continued n Page 70

Read and Understood By

J. Cabralda

Signed

Date

Quis-

Signed

Dat

- selected 3X clones #1, 12 and 4X clones #2, 10 to work with further
- electroporated each of these into MC4100(+), using 1 μ l plasmid DNA and 1 μ l pGP1-2, as per SOP
- plated 10 μ l of each \rightarrow incubated 30°C O/N
(Mac + Kan + Amp)
- very few colonies!
- only 2 plates (#10 & #12) have any colonies - and only a few colonies each!
- luckily, this includes a 3X and 4X plate so I can continue with the experiment

PREPARATION OF OVERNIGHT CULTURES

Overnights

- prepared two culture tubes, each with 3mL TB + 100 μ g/mL Ampicillin + 30 μ g/mL Kanamycin
- inoculated one tube with #10 and the other with #12 (using toothpicks)
- grew both tubes O/N @ 30°C, 300 rpm.

VERIFICATION OF PRESENCE OF PLASMIDS

Plasmid Preps

- removed 1.8 mL of each overnight culture to a fresh microfuge tube
- did plasmid preps as per SOP
- stopped after isopropanol ppt \rightarrow resuspended in 100 μ l TE + RNase

Gels

- prepared a 1% agarose gel

Continued on Page

72

Read and Understood By

J. Cabralda

Signed

Date


Opini

Signed

Date

- loaded 2 μ l each of : pGP1-2, #12, #10, #12 w/pGP, #10 w/pGP
- ran gel 60' @ 60V.

Figure 3 Plasmid preps of putative transformants. Both are positive for pGP1-2.

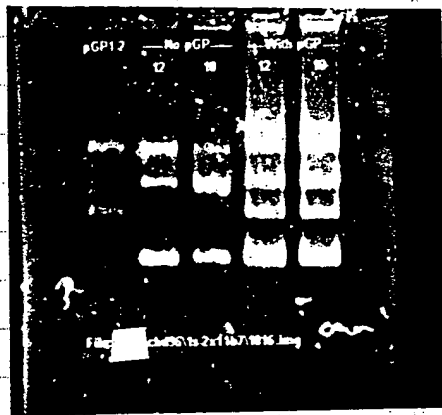

INDUCTION OF EXPRESSION

Harvesting of Preinduction Cultures

- from each O/N culture, removed 100 μ l to a fresh microfuge tube
- @ 13,000 rpm, 1' \rightarrow removed S/N
- resuspended each pellet in 100 μ l SDS reducing buffer
- froze @ -20°C until needed.

Induction of Culture

- removed 1ml of each O/N culture to a fresh culture tube containing 1ml fresh TB
- incubated 3h @ 42°C, with shaking (300 rpm)

CHECKING FOR FUSION PROTEIN EXPRESSION

Continued on Page

73

Read and Understood By

J. Cabralda

Signed

Date

[Signature]

Signed

Date

Harvesting of Induction Culture

- removed 200 μ l of each culture to microfuge tubes
- ϕ 11' \rightarrow removed SYN
- resuspended each pellet in 100 μ l SDS reducing buffer


Preparation of Samples

- removed preinduced samples from freezer
- boiled all samples 10'
- spun down to concentrate

Gel

- prepared a 15% / 4% stacking gel
- loaded 10 μ l each sample, along with 5 μ l low MW standards.
- ran gel @ 100V for 2h then @ 200V for 10'

Figure 4 Expression of CBD96-1S-3x11B7 and CBD96-2S-4x11B7 fusion proteins.

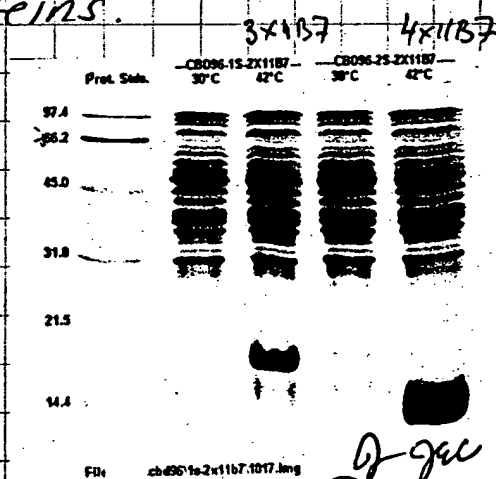

- expression of 1S-2x11B7 is excellent!
- but 2S-2x11B7 is far too small - perhaps this clone has a frameshift mutation close to CBD96?
- or maybe it has no insert after all

Continued on Pag

75

Read and Understood By

J. Cabralda

Signed

Date

J. Cabralda

Sign d

Dat